



BLUE
PORTABLE ULTRASOUND SCANNER
FOR COMPLETE EQUINE DIAGNOSTICS


BLUE

PORTABLE ULTRASOUND SCANNER

A premium class mobile ultrasound machine for quick, precise and comprehensive health diagnosis in horses. Intuitive menu, crystal clear touch screen and user friendly interface create a great tool for high quality imaging both in clinic and field environment.

BLUE FEATURES

- Responsive and crystal clear touchscreen for high quality imaging. Facilitates differentiation of examined structures beyond reproductive tract (lungs, udder, liver, etc.).
- Short boot-up time and active fields for faster set up.
- User friendly interface for better workflow.
- Compatible with tall stand and tendon stand for higher mobility and flexibility of use.
- Remote training and after-sales service for peace of mind.

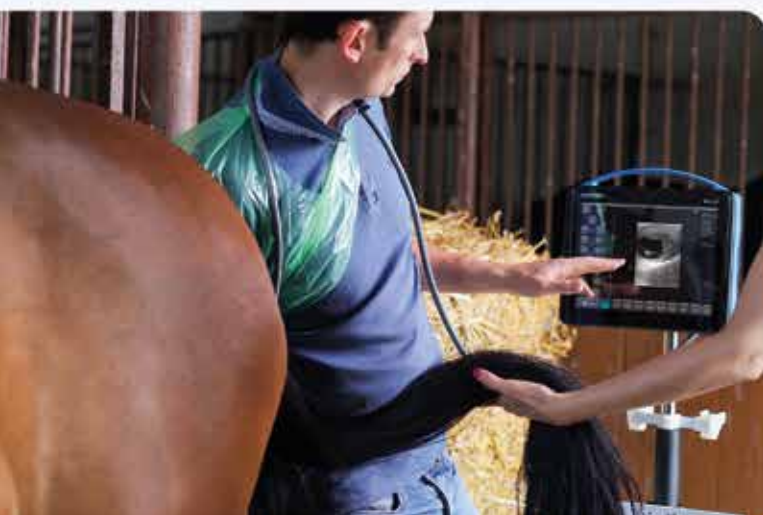



FEATURES OF BLUE



12" LED LCD touch screen.



Weighs only 4 kg.



Internal battery lasting 2.5 hours of continuous work and the option of external power supply.



Sensitive color Doppler for easy foetus heartbeat detection.






Tendon stand




Musculoskeletal system



Pregnancy

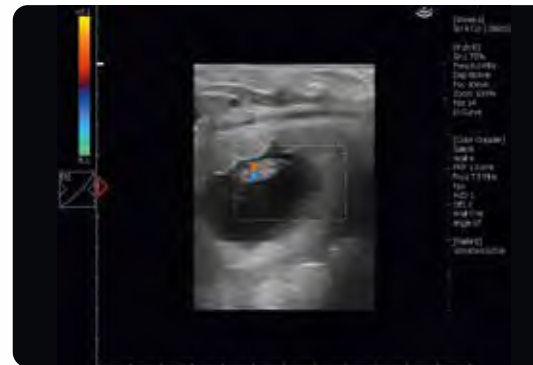

Pregnancy (Doppler image)

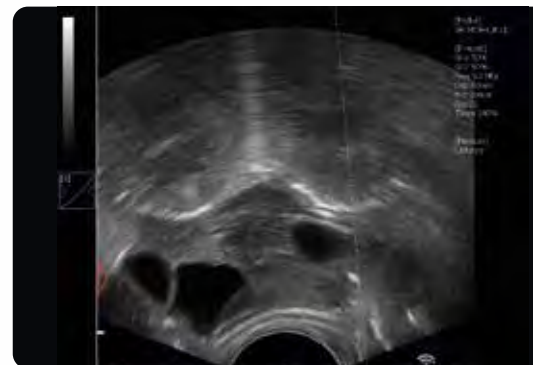

Ovary - OPU

TECHNOLOGIES OF THE FUTURE:



for linear probes allows a wider field of view,



grayscale curve perfectly matched to the perception of the human eye,



contrast, sharpness and tissue imaging enhancement system.



PROBES AND ACCESSORIES

Linear abdominal

frequency:
6 - 14 MHz

scanning depth:
40 - 100 mm

footprint:
40 mm



The OPU system

Fast to assembly and disassembly.
Made from lightweight and easy to clean material.



Convex abdominal

frequency:
2 - 8 MHz

scanning depth:
60 - 300 mm

footprint:
50 mm



WiFi antenna (USB)

Enables connecting Dрамиński BLUE ultrasound scanner with Internet.



Convex abdominal

frequency:
2 - 5 MHz

scanning depth:
60 - 350 mm

footprint:
60 mm



Gel Standoffs

Enable examination of tendons with the use of a linear probe.



Linear rectal

frequency:
4 - 9 MHz

scanning depth:
40 - 150 mm

footprint:
60 mm



Sun protection shield

Facilitates work in strong lighting.





SEE MORE

DRAMIŃSKI S.A.

Wiktora Steffena 21, 11-036 Sząbruk, Poland

 +48 89 675 26 00  ultrasound@draminski.com

www.draminski.com